# ABclonal www.abclonal.com

# Calponin Rabbit mAb

Catalog No.: A10822 Recombinant

# **Basic Information**

#### **Observed MW**

33kDa

#### **Calculated MW**

33kDa

#### **Category**

SMab Recombinant Monoclonal Antibody

## **Applications**

WB,IHC-P,IF/ICC,ELISA

### **Cross-Reactivity**

Human, Mouse, Rat, Dog, Pig

# **Background**

Predicted to enable actin binding activity. Involved in negative regulation of vascular associated smooth muscle cell proliferation. Located in cytoskeleton.

# **Recommended Dilutions**

WB	1:500 - 1:2000
IHC-P	1:50 - 1:200
IF/ICC	1:50 - 1:200

# **Immunogen Information**

Gene ID	<b>Swiss Prot</b>
1264	P51911

#### **Immunogen**

Recombinant protein of human Calponin

#### **Synonyms**

SMCC; Sm-Calp; HEL-S-14; Calponin

### Contact


# **Product Information**

Source	Isotype	Purification
Rabbit	IqG	Affinity purification

# Storage

Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

# Validation Data

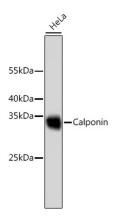

Western blot analysis of lysates from HeLa cells, using Calponin Rabbit mAb (A10822). Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins:  $25\mu g$  per lane.

Blocking buffer: 3% nonfat dry milk in TBST.