

Recombinant Human Follicle-Stimulating Hormone α/β Dimer Protein

Catalog No.: RP00702 Recombinant

Sequence Information

Species Gene ID Swiss Prot Human P01215&P0 1225

Tags No tag

Synonyms

Follicle-stimulating hormone; FSH; FSH alpha/beta

Product Information

Purification Source HEK293 cells > 95% by SDS-

PAGE.

Endotoxin

< 1 EU/µg of the protein by LAL method.

Formulation

Lyophilized from a 0.2 µm filtered solution of PBS, pH7.4.Contact us for customized product form or formulation.

Reconstitution

Reconstitute to a concentration of 0.1-0.5 mg/mL in sterile distilled water.

Background

Human Follicle-stimulating hormone (FSH) is a member of glycoprotein hormones subunit beta family, which also includes LH, chorionic gonadotropin (CG) and thyroid-stimulating hormone (TSH). FSH and itsfamilymembers are heterodimers consisting of non-covalently linked α - and β -subunits. They share an identical α subunit, and β -subunits vary. FSH has a unique β -subunit (FSH β), which confers its specific biologic activityandis responsible for interaction with the FSH-receptor which belongs to a subfamily of GPCRs calledleucine-rich-repeatcontaining GPCRs (LGRs). FSH is secreted from the pituitary gland and regulatesreproduction inmammals. FSH stimulates sertoli cell proliferation in testes and supports spermatogenesis inmales, andinduces the maturation of ovarian follicles in females.

Basic Information

Description

Recombinant Human Follicle-Stimulating Hormone α/β Dimer Protein is produced by Human Cells expression system. The target protein is expressed with sequence (Ala25-Ser116&Asn19-Glu129) of human Follicle-Stimulating Hormone α/β Dimer (Accession #P01215&P01225).

Bio-Activity

Storage

Store the lyophilized protein at -20°C to -80 °C for long term.
 After reconstitution, the protein solution is stable at -20 °C for 3 months, at 2-8 °C for up to 1 week.


Avoid repeated freeze/thaw cycles.

Contact



www.abclonal.com

Validation Data

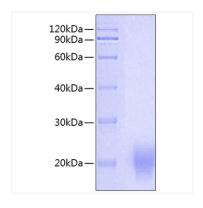

Recombinant Human Follicle-Stimulating Hormone α/β Dimer Protein was determined by SDS-PAGE under reducing conditions with Coomassie Blue.