# ABclonal www.abclonal.com

## PE Rabbit anti-Human CD85j/ILT2 mAb

Catalog No.: A24250

## **Basic Information**

#### **Observed MW**

## **Calculated MW**

49kDa/70kDa/71kDa

## Category

SMab Recombinant Monoclonal Antibody

#### **Applications**

FC

## **Cross-Reactivity**

Human

#### CloneNo number

ARC60700-PE

## Conjugate

PE. Ex:565nm. Em:574nm.

## Background

This gene is a member of the leukocyte immunoglobulin-like receptor (LIR) family, which is found in a gene cluster at chromosomal region 19q13.4. The encoded protein belongs to the subfamily B class of LIR receptors which contain two or four extracellular immunoglobulin domains, a transmembrane domain, and two to four cytoplasmic immunoreceptor tyrosine-based inhibitory motifs (ITIMs). The receptor is expressed on immune cells where it binds to MHC class I molecules on antigen-presenting cells and transduces a negative signal that inhibits stimulation of an immune response. It is thought to control inflammatory responses and cytotoxicity to help focus the immune response and limit autoreactivity. Multiple transcript variants encoding different isoforms have been found for this gene.

## **Recommended Dilutions**

FC

5 μl per 10^6 cells in 100 μl volume

## **Immunogen Information**

**Gene ID** 10859

Swiss Prot Q8NHL6

#### **Immunogen**

Recombinant fusion protein containing a sequence corresponding to amino acids 24-458 of human CD85j/ILT2 (NP\_006660.4).

### **Synonyms**

LILRB1; CD85J; ILT-2; ILT2; LIR-1; LIR1; MIR-7; MIR7; PIR-B; PIRB; leukocyte immunoglobulin like receptor B1

#### **Contact**

•

www.abclonal.com

## **Product Information**


SourceIsotypePurificationRabbitIgGAffinity purification

## Storage


Store at 2-8°C. Avoid freeze.

Buffer: PBS with 0.03% proclin300,0.2% BSA,pH7.3.

## **Validation Data**

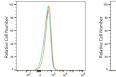



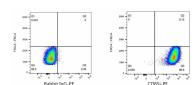

Flow cytometry: 1X10^6 K-562 cells(negative control,left) and Daudi cells(right) were surface-stained with PE Rabbit anti-Human CD85j/ILT2 mAb(A24250,5 µl/Test,orange line) or PE Rabbit IgG isotype control (5 µl/Test,blue line).Non-fluorescently stained cells were used as blank control (red line).

Flow cytometry:1X10^6 Daudi cells were surface-stained with PE Rabbit IgG isotype control (5  $\mu$ I/Test,left) or PE Rabbit anti-Human CD85j/ILT2 mAb(A24250,5  $\mu$ I/Test,right).