

AS002

Leader in Biomolecular Solutions for Life Science



HRP-conjugated Goat anti-Human IgG (H+L)

Catalog No.: AS002 **20 Publications**

Basic Information

Observed MW

25kDa,55kDa

Calculated MW

Category

Antibody

Applications

WB

Cross-Reactivity

Conjugate

HRP

Background

Secondary antibodies are affinity-purified antibodies which will work with target-specific primary antibody in the detection, sorting or purification of its specified target. Secondary antibodies offer increased versatility enabling users to use many detection systems (e.g. HRP, AP, fluorescence). They can also provide greater sensitivity through signal amplification as multiple secondary antibodies. Most commonly, secondary antibodies are generated by immunizing the host animal (different from host species of primary antibody) with a pooled population of normal immunoglobulins from the host species of primary antibody and can be further purified and modified (i.e. antibody fragmentation, label conjugation, etc.) to ensure well-characterized specificity to corresponding normal immunoglobulins.

Recommended Dilutions

WB 1:1000 - 1:10000

Immunogen Information

Gene ID Swiss Prot

Immunogen

Human IgG

Synonyms

Contact

 www.abclonal.com

Product Information

Source

Goat

Isotype

Horseradish peroxidase
conjugated IgG


Purification

Affinity purification


Storage

Store at -20°C. Avoid freeze / thaw cycles.
Buffer: PBS with 0.75% BSA,50% glycerol,pH7.3.

Validation Data

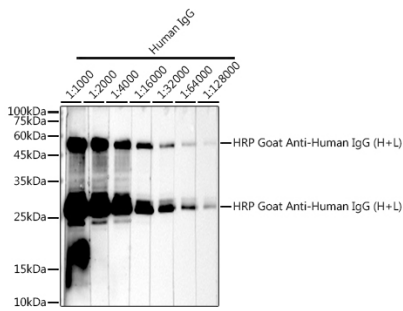

Western blot analysis of lysates from Human IgG, using HRP Goat Anti-Human IgG (H+L) (AS002) at 1:2000 dilution.
Lysates/proteins: 25µg per lane.
Blocking buffer: 3% nonfat dry milk in TBST.
Detection: ECL Basic Kit (RM00020).
Exposure time: 90s.

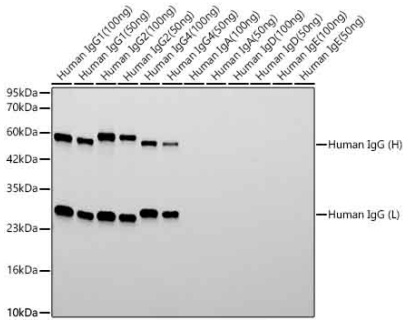

Western blot analysis of lysates from Human IgG using HRP Goat Anti-Human IgG (H+L) (AS002) at 1:10000 dilution.
Lysates/proteins: 50ng - 100ng per lane.
Blocking buffer: 3% nonfat dry milk in TBST.
Detection: ECL Basic Kit (RM00020).
Exposure time: 20s.