

Catalog No.: A10536 2 Publications



Basic Information

Observed MW 75kDa

Calculated MW 75kDa

Category Polyclonal Antibody

Applications WB,IHC-P,IF/ICC,ELISA

Cross-Reactivity Human, Mouse, Rat

Background

DVL1, the human homolog of the Drosophila dishevelled gene (dsh) encodes a cytoplasmic phosphoprotein that regulates cell proliferation, acting as a transducer molecule for developmental processes, including segmentation and neuroblast specification. DVL1 is a candidate gene for neuroblastomatous transformation. The Schwartz-Jampel syndrome and Charcot-Marie-Tooth disease type 2A have been mapped to the same region as DVL1. The phenotypes of these diseases may be consistent with defects which might be expected from aberrant expression of a DVL gene during development.

Recommended Dilutions

WB	1:500 - 1:1000
IHC-P	1:50 - 1:200
IF/ICC	1:50 - 1:200
ELISA	Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID	Swiss Prot	
1855	O14640	
Immunogen Recombinant fusion protein containing a sequence corresponding to amino acids 510-640 of human DVL1 (NP_004412.2).		
Synonyms DVL; DRS2; DVL1L1; DVL1P1; DVL1		
Product Information		

Contact

• www.abclonal.com


Source Rabbit

Isotype IgG **Purification** Affinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.

Validation Data

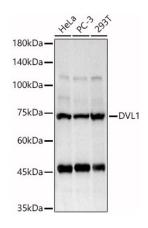

Western blot analysis of various lysates, using DVL1 Rabbit pAb (A10536) at 1:1000 dilution.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 20s.

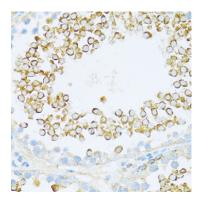


Immunohistochemistry analysis of paraffin-embedded Rat testis using DVL1 Rabbit pAb (A10536) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

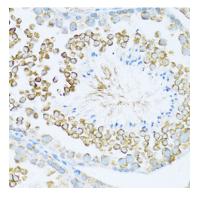

Immunohistochemistry analysis of paraffin-embedded Mouse testis using DVL1 Rabbit pAb (A10536) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

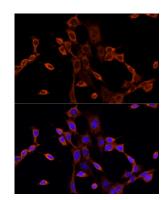

Immunofluorescence analysis of NIH-3T3 cells using DVL1 Rabbit pAb (A10536) at dilution of 1:100 (40x lens). Secondary antibody: Cy3conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

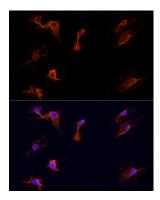

Immunofluorescence analysis of C6 cells using DVL1 Rabbit pAb (A10536) at dilution of 1:100 (40x lens). Secondary antibody: Cy3conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.