ABclonal www.abclonal.com

HCoV-229E Spike S1 Rabbit pAb

Catalog No.: A20393

Basic Information

Observed MW

160kDa

Calculated MW

129kDa

Category

Polyclonal Antibody

Applications

WB, ELISA

Cross-Reactivity

HCoV-229E

Background

S1 region attaches the virion to the cell membrane by interacting with host ANPEP/aminopeptidase N, initiating the infection. Binding to the receptor probably induces conformational changes in the S glycoprotein unmasking the fusion peptide of S2 region and activating membranes fusion. S2 region belongs to the class I viral fusion protein. Under the current model, the protein has at least 3 conformational states: prefusion native state, pre-hairpin intermediate state, and post-fusion hairpin state. During viral and target cell membrane fusion, the coiled coil regions (heptad repeats regions assume a trimer-of-hairpins structure, positioning the fusion peptide in close proximity to the C-terminal region of the ectodomain. The formation of this structure appears to drive apposition and subsequent fusion of viral and target cell membranes.

Recommended Dilutions

WB 1:500 - 1:1000

ELISA

Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific

assay requirements.

Immunogen Information

Gene ID918758

Swiss Prot
P15423

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 1074-1173 of coronavirus Spike S1 (NP_073551.1).

Synonyms

Contact

www.abclonal.com


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.

Validation Data

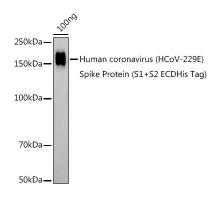

Western blot analysis of lysates from Human coronavirus (HCoV-229E) Spike Protein (S1+S2 ECDHis Tag), using HCoV-229E Spike S1 Rabbit pAb (A20393) at 1:1000 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 180s.