


6x DNA Loading Dye

#Cat: NB-03-0169 Size: 5x1ml

Contents

10 mM Tris-HCL (pH 7.6) 0.03 % Bromophenol Blue 0.03% Xylene Cyanol 60 % Glycerol 60 mM EDTA

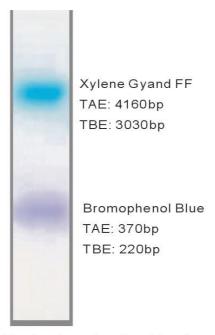

Electrophoresis of tracking dyes in 6x Gel Loading Dye

Description

The Neobiotech 6x DNA Loading Dye is a pre-mixed loading buffer with two tracking dye for agarose and non-denaturing poylacrylamide gel electrophoresis. EDTA is included to chelate magnesium (up to 10 mM) in enzymatic reactions, thereby stopping the reaction. Bromophenol blue and xylene cyanol are the standard tracking dye for electrophoresis.

Applications

Preparation of DNA ladders, markers and samples for loading on agarose or polyacrylamide gels.

Features

Two-color tracking of DNA migration during electrophoresis. No

DNA masking during gel exposure to UV light.

EDTA binds divalent metal ions and inhibits metal dependent nucleases.



Storage

1 year at room temperature, for long-time preservation, please storage at -20°C.

Quality Assurance Statement

Gel Loading Dye(6X) is assayed for contaminating endonuclease, exonuclease and RNase Activity