

RAB11FIP1 antibody

Product Information

Catalog No.: FNab06990

Size: 100μg Form: liquid

Purification: Immunogen affinity purified

Purity: ≥95% as determined by SDS-PAGE

Host: Rabbit

Clonality: polyclonal

Clone ID: None IsoType: IgG

Storage: PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12

months(Avoid repeated freeze / thaw cycles.)

Background

A Rab11 effector protein involved in the endosomal recycling process. Also involved in controlling membrane trafficking along the phagocytic pathway and in phagocytosis.

Immunogen information

Immunogen: RAB11 family interacting protein 1(class I)

Synonyms: Rab11 family-interacting protein 1 (Rab11-FIP1)|Rab-coupling

protein|RAB11FIP1|RCP

Observed MW: 71 kDa
Uniprot ID: Q6WKZ4

Application


Reactivity: Human, Mouse, Rat, Monkey

Tested Application: ELISA, WB, IHC, IF


Recommended dilution: WB: 1:500-1:2000; IHC: 1:20-1:200; IF: 1:20-1:200

Image:



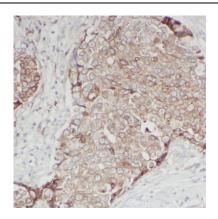

Immunohistochemistry of paraffin-embedded human breast cancer tissue slide using FNab06990(RAB11FIP1 Antibody) at dilution of 1:50

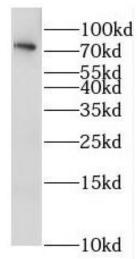

HepG2 cells were subjected to SDS PAGE followed by western blot with FNab06990(RAB11FIP1 Antibody) at dilution of 1:600