
phospho-TDP43 (403/404) antibody

Product Information

Catalog No.:	FNab06400
Size:	100µg
Form:	liquid
Purification:	Protein A+G purification
Purity:	≥95% as determined by SDS-PAGE
Host:	Mouse
Clonality:	monoclonal
Clone ID:	0F8
IsoType:	IgG1
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

Background


This antibody is a mouse monoclonal antibody recognizing TDP-43 only when phosphorylated at 403/404. This antibody detects phosphorylated forms of TDP-43 as well as the cleavage product of 20-30 kDa.

Immunogen information

Immunogen:	TAR DNA binding protein
Synonyms:	TAR DNA-binding protein 43 (TDP-43) TARDBP TDP43
Observed MW:	25 kDa
Uniprot ID :	Q13148

Application

Reactivity:	Human, Mouse
Tested Application:	ELISA, WB
Recommended dilution:WB:	1:1000-1:6000
Image:	

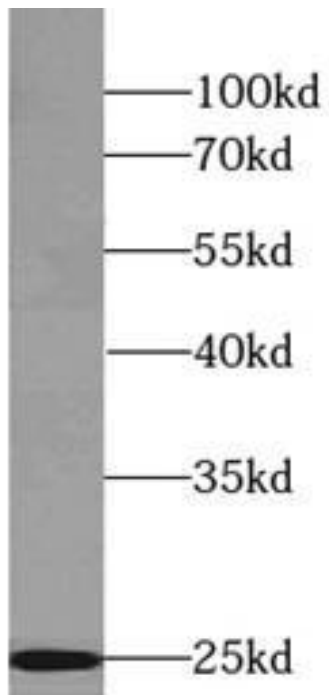

K-562 cells were subjected to SDS PAGE followed by western blot with FNab06400(phospho(403/404)-TDP43 Antibody) at dilution of 1:3000