

Human Thioredoxin Reductase 1/TRXR1 Rabbit mAb

Catalog No.: A4725 Recombinant 1 Publications

Basic Information

Observed MW

60kDa, 71kDa

Calculated MW

71kDa

Category

SMab Recombinant Monoclonal Antibody

Applications

WB, ELISA

Cross-Reactivity

Human, Mouse, Rat

CloneNo number

ARC1106

Background

The protein encoded by this gene belongs to the pyridine nucleotide-disulfide oxidoreductase family, and is a member of the thioredoxin (Trx) system. Three thioredoxin reductase (TrxR) isozymes are found in mammals. TrxRs are selenocysteine-containing flavoenzymes, which reduce thioredoxins, as well as other substrates, and play a key role in redox homoeostasis. This gene encodes an ubiquitously expressed, cytosolic form of TrxR, which functions as a homodimer containing FAD, and selenocysteine (Sec) at the active site. Sec is encoded by UGA codon that normally signals translation termination. The 3' UTRs of selenoprotein mRNAs contain a conserved stem-loop structure, the Sec insertion sequence (SECIS) element, which is necessary for the recognition of UGA as a Sec codon rather than as a stop signal. Alternative splicing, primarily at the 5' end, results in transcript variants encoding same or different isoforms, including a glutaredoxin-containing isoform that is predominantly expressed in testis.

Recommended Dilutions

WB 1:1000 - 1:4000

ELISA

Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID7296
Swiss Prot
Q16881

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 1-649 of human Thioredoxin reductase 1 (TXNRD1) (Q16881).

Synonyms

TR; TR1; TXNR; TRXR1; GRIM-12; Thioredoxin reductase 1 (TXNRD1)

Contact

www.abclonal.com

Product Information


SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

Validation Data

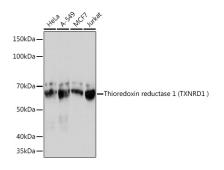

Western blot analysis of various lysates using Thioredoxin reductase 1 (Thioredoxin reductase 1 (TXNRD1)) Rabbit mAb (A4725) at 1:1000 dilution.

Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000


dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 10s.

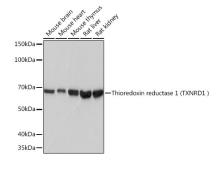

Western blot analysis of various lysates using Thioredoxin reductase 1 (Thioredoxin reductase 1 (TXNRD1)) Rabbit mAb (A4725) at 1:1000 dilution.

Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000

dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 30s.