ABclonal www.abclonal.com

YAP1 Rabbit pAb

Catalog No.: A11918

Basic Information

Observed MW

65kDa

Calculated MW

54kDa

Category

Polyclonal Antibody

Applications

WB,IF/ICC,IP,ELISA

Cross-Reactivity

Human, Mouse, Zebrafish

Background

This gene encodes a downstream nuclear effector of the Hippo signaling pathway which is involved in development, growth, repair, and homeostasis. This gene is known to play a role in the development and progression of multiple cancers as a transcriptional regulator of this signaling pathway and may function as a potential target for cancer treatment. Alternative splicing results in multiple transcript variants encoding different isoforms.

Recommended Dilutions

WB	1:500 - 1:2000
IF/ICC	1:50 - 1:100
ΤD	1·20 - 1·50

Immunogen Information

Gene ID	Swiss Prot
10413	P46937

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 350-450 of human YAP1 (NP_001123617.1).

Synonyms

YAP; YKI; COB1; YAP2; YAP-1; YAP65; YAP1

Contact

\odot	www.abclonal.com

Product Information


SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

Validation Data

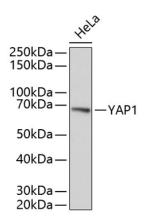

Western blot analysis of lysates from HeLa cells, using YAP1 Rabbit pAb (A11918). Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.