

A2202

Leader in Biomolecular Solutions for Life Science



TOLLIP Rabbit pAb

Catalog No.: A2202 **4 Publications**

Basic Information

Observed MW

30kDa

Calculated MW

30kDa

Category

Polyclonal Antibody

Applications

WB,IHC-P,IF/ICC,ELISA

Cross-Reactivity

Human,Mouse,Rat

Background

This gene encodes a ubiquitin-binding protein that interacts with several Toll-like receptor (TLR) signaling cascade components. The encoded protein regulates inflammatory signaling and is involved in interleukin-1 receptor trafficking and in the turnover of IL1R-associated kinase. Several transcript variants encoding different isoforms have been found for this gene.

Recommended Dilutions

WB 1:500 - 1:1000

IHC-P 1:50 - 1:200

IF/ICC 1:50 - 1:200

ELISA Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

Immunogen Information

Gene ID

54472

Swiss Prot

Q9H0E2

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 1-274 of human TOLLIP (NP_061882.2).

Synonyms

IL-1RAcPIP; TOLLIP

Contact



www.abclonal.com

Product Information

Source

Rabbit

Isotype

IgG

Purification


Affinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.09% Sodium azide,50% glycerol,pH7.3.

Validation Data

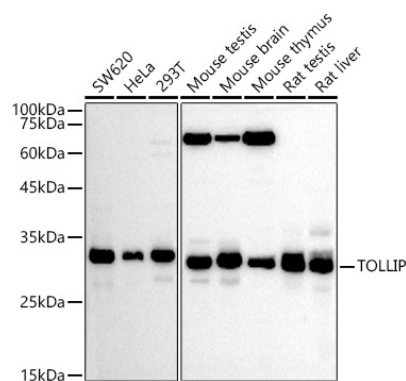

Western blot analysis of various lysates using TOLLIP Rabbit pAb (A2202) at 1:500 dilution.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.


Detection: ECL Basic Kit (RM00020).

Exposure time: 30s.

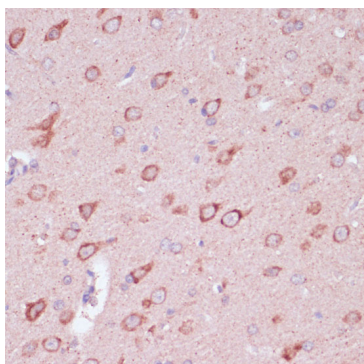

Immunohistochemistry analysis of paraffin-embedded Rat brain using TOLLIP Rabbit pAb (A2202) at dilution of 1:100 (40x lens).

Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

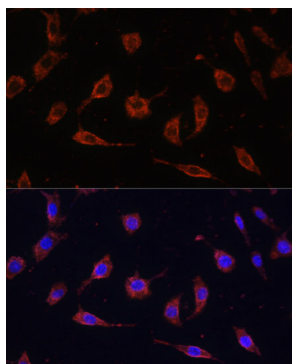

Immunofluorescence analysis of L929 cells using TOLLIP Rabbit pAb (A2202) at dilution of 1:100.

Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.