Leader in Biomolecular Solutions for Life Science

# **AMSH Rabbit mAb**

Catalog No.: A20889 Recombinant



### **Basic Information**

### **Observed MW**

50kDa

### **Calculated MW**

48kDa

#### Category

SMab Recombinant Monoclonal Antibody

### **Applications**

WB, ELISA

### **Cross-Reactivity**

Human, Monkey

### CloneNo number

ARC2840

# **Background**

Cytokine-mediated signal transduction in the JAK-STAT cascade requires the involvement of adaptor molecules. One such signal-transducing adaptor molecule contains an SH3 domain that is required for induction of MYC and cell growth. The protein encoded by this gene binds to the SH3 domain of the signal-transducing adaptor molecule, and plays a critical role in cytokine-mediated signaling for MYC induction and cell cycle progression. Multiple alternatively spliced transcript variants encoding the same protein isoform have been found for this gene.

# **Recommended Dilutions**

WB

1:500 - 1:1000

# **Immunogen Information**

**Gene ID** 10617

Swiss Prot 095630

#### **Immunogen**

A synthetic peptide corresponding to a sequence within amino acids 100-200 of human AMSH (095630).

### **Synonyms**

AMSH; MICCAP

### **Contact**

•

www.abclonal.com

# **Product Information**


**Source** Rabbit **Isotype** IgG **Purification**Affinity purification

# Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

# **Validation Data**

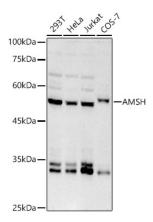

Western blot analysis of extracts of various cell lines, using AMSH antibody (A20889) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane. Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 180s.