


## Antibodies - Chicken Ig Why?

Immunoglobulin Y (IgY) is the dominant immunoglobulin found in avian, amphibian and reptilian species. In 1893, Klemperer[1] first described the concept of passive immunity by reporting the presence and transference of tetanus toxin immunity (antibodies) from hen to chick via egg. He noted that antibodies were present in both the blood of the chicken and the yolk of the egg that it laid. However, it wasn't until 1969 that Leslie and Clem[2] coined the term IgY to describe poultry antibodies, including those found in the yolk of the eggs. Until 1969, IgY was known and reported as IgG, but the structural and antigenic differences were shown to be great enough to earn the protein the new Immunoglobulin classification IgY (Figure 1).

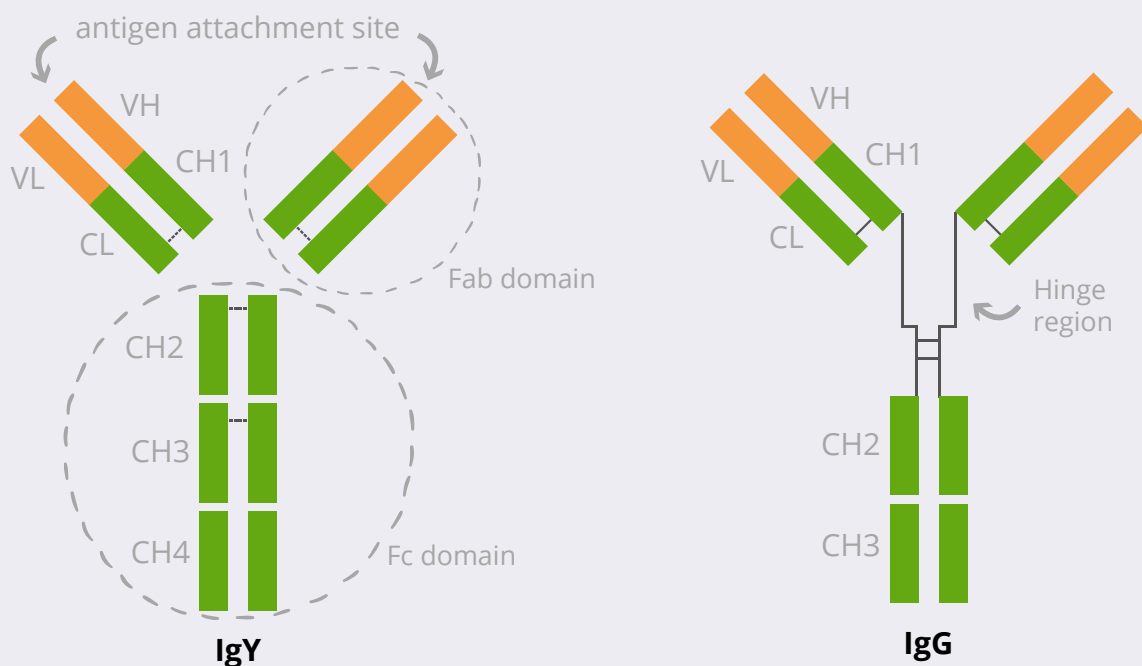

(Figure 1) Avian IgY structure, compared to mammalian IgG. Both contain two heavy and two light chains, which consist of a variable domain (VH and VL) and four constant domains (CH1, CH2, CH3 and CH4). IgG has a longer hinge region making it more flexible than avian IgY

### IgG and IgY are similar on several levels:

- Both have the characteristic antibody Y shape
- Both contain 2 heavy and 2 light chains
- Both are divalent
- Both play a similar biological role, being the major immunoglobulins to provide defense against infectious agents, appearing in high concentrations after initial synthesis of the higher molecular weight antibody (IgM)

### IgG and IgY are different on several levels:

- Little or no immunological cross-reactivity between IgY and mammalian IgG
- IgY has a higher molecular weight due to an extra heavy chain constant domain
- IgY lacks a well-defined hinge region
- IgY heavy chains lack Fc domains so they don't fix complement or bind protein A or protein G
- IgY can facilitate anaphylactic reactions, which can only be done by IgE in mammals




IgY also shares some similarities with mammalian IgE, including similar intrachain disulphide bonding in their extra heavy chain domain. This led to suggestions that IgY may well be an ancestral molecule to both mammalian IgG and IgE, which have since been confirmed by genetic and structural studies[3][4][5].

**Chicken  
raised  
antibodies**



**IgY**

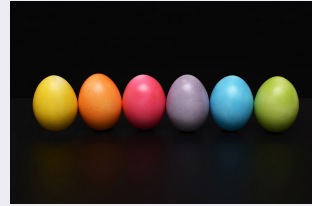

## Why would you want to use IgY's in your research?

The very fact that the molecules are so similar in both their biological role and their structure, yet different enough not to cross-react offers several advantages to a research scientist. In most experimental applications, including Western blotting, immunohistochemistry, immunocytochemistry, ELISA and functional blocking experiments, both IgG and IgY are equivalent, meaning IgY can generally be used as a direct alternative to IgG[5].

## Advantages of IgY

- **Higher avidity** - Chickens are not mammals so their reaction to a mammalian antigen is enhanced leading to high-avidity antibodies. The response is particularly strong to mammalian antigens and especially those which are highly conserved.
- **Humane** - What could be more humane? The chicken lays the egg, so there is no need to bleed or otherwise conduct procedures on the chicken; you just collect the eggs. Plus, Chickens are cheaper to keep than rabbits![6]
- **High yield** - one chicken can produce large quantities of IgY against low quantities of antigen. Some can produce 2.5 - 3 g of IgY per month, which is 10-20 times the amount of a rabbit[6]
- **Quick** - Compared to rabbits, you can obtain high-titre antibody from eggs as early as day 25.
- **Lower background** - IgY does not bind to IgG Fc receptors meaning you get less false positive staining. They also do not activate mammalian complement systems.
- **Cross-reactivity** - As there is none to mammalian species, they make an excellent addition to multiple labelling or multi-plex experiments.
- **Convenient** - IgY is packed neatly in eggs, which can be stored for long periods and the IgY purified to the desired titre/avidity months later.
- IgY can still be digested by papain to produce a divalent Fab fragment.
- IgY can be conjugated to enzymes or fluorophores, biotinylated and gold-labeled by standard procedures, offering you the full range of functionality you would expect from an antibody tool.

## Helping your Research

So, why not add these fantastic tools to your kit when you next need to visualize a protein? Biorbyt has over 400 Chicken raised antibodies to offer, so check out our website to see if we have added your target to our portfolio.

## References

1. Klemperer, F. Ueber natürliche Immunität und ihre Verwerthung für die Immunisierungstherapie. Arch. Für Exp. Pathol. Pharmakol. 31, 356-382 (1893).
2. Leslie, G. A. & Clem, L. W. Phylogen of immunoglobulin structure and function. 3. Immunoglobulins of the chicken. J. Exp. Med. 130, 1337-1352 (1969).
3. Zhang, X., Calvert, R. A., Sutton, B. J. & Doré, K. A. IgY: a key isotype in antibody evolution: IgY antibody. Biol. Rev. 92, 2144-2156 (2017).
4. Taylor, A. I., Fabiane, S. M., Sutton, B. J. & Calvert, R. A. The Crystal Structure of an Avian IgY-Fc Fragment Reveals Conservation with both Mammalian IgG and IgE † ‡. Biochemistry 48, 558-562 (2009).
5. Lee, W., Syed Atif, A., Tan, S. C. & Leow, C. H. Insights into the chicken IgY with emphasis on the generation and applications of chicken recombinant monoclonal antibodies. J. Immunol. Methods 447, 71-85 (2017).
6. Amro, W. A., Al-Qaisi, W. & Al-Razem, F. Production and purification of IgY antibodies from chicken egg yolk. J. Genet. Eng. Biotechnol. 16, 99-103 (2018).

**Distributed by:**

**CliniSciences Group**

# CliniSciences Group

## Austria

Company: CliniSciences GmbH  
Address: Sternwartestrasse 76, A-1180  
Wien - Austria  
Telephone: +43 720 115 580  
Fax: +43 720 115 577  
Email: [oesterreich@clinisciences.com](mailto:oesterreich@clinisciences.com)  
Web: <https://www.clinisciences.com>



## Belgium

Company: CliniSciences S.R.L  
Address: Avenue Stalingrad 52, 1000  
Brussels - Belgium  
Telephone: +32 2 31 50 800  
Fax: +32 2 31 50 801  
Email: [belgium@clinisciences.com](mailto:belgium@clinisciences.com)  
Web: <https://www.clinisciences.com>



## Denmark

Company: CliniSciences ApS  
Address: Oesterbrogade 226, st. 1,  
Copenhagen, 2100 - Denmark  
Telephone: +45 89 888 349  
Fax: +45 89 884 064  
Email: [danmark@clinisciences.com](mailto:danmark@clinisciences.com)  
Web: <https://www.clinisciences.com>



## Finland

Company: CliniSciences ApS  
Address: Oesterbrogade 226, st. 1,  
Copenhagen, 2100 - Denmark  
Telephone: +45 89 888 349  
Fax: +45 89 884 064  
Email: [suomi@clinisciences.com](mailto:suomi@clinisciences.com)  
Web: <https://www.clinisciences.com>



## France

Company: CliniSciences S.A.S  
Address: 74 Rue des Suisses, 92000  
Nanterre- France  
Telephone: +33 9 77 40 09 09  
Fax: +33 9 77 40 10 11  
Email: [info@clinisciences.com](mailto:info@clinisciences.com)  
Web: <https://www.clinisciences.com>



## Germany

Company: Biotrend Chemikalien GmbH  
Address: Wilhelm-Mauser-Str. 41-43,  
50827 Köln - Germany  
Telephone: +49 221 9498 320  
Fax: +49 221 9498 325  
Email: [info@biotrend.com](mailto:info@biotrend.com)  
Web: <https://www.biotrend.com>



## Iceland

Company: CliniSciences ApS  
Address: Oesterbrogade 226, st. 1,  
Copenhagen, 2100 - Denmark  
Telephone: +45 89 888 349  
Fax: +45 89 884 064  
Email: [island@clinisciences.com](mailto:island@clinisciences.com)  
Web: <https://www.clinisciences.com>



## Ireland

Company: CliniSciences Limited  
Address: Ground Floor, 71 lower Baggot street  
Dublin D02 P593 - Ireland  
Telephone: +353 1 6971 146  
Fax: +353 1 6971 147  
Email: [ireland@clinisciences.com](mailto:ireland@clinisciences.com)  
Web: <https://www.clinisciences.com>



## Italy

Company: CliniSciences S.r.l  
Address: Via Maremmana inferiore 378  
Roma 00012 Guidonia Montecelio - Italy  
Telephone: +39 06 94 80 56 71  
Fax: +39 06 94 80 00 21  
Email: [italia@clinisciences.com](mailto:italia@clinisciences.com)  
Web: <https://www.clinisciences.com>



## Netherlands

Company: CliniSciences B.V.  
Address: Kraaijenhoffstraat 137A,  
1018RG Amsterdam, Netherlands  
Telephone: +31 85 2082 351  
Fax: +31 85 2082 353  
Email: [nederland@clinisciences.com](mailto:nederland@clinisciences.com)  
Web: <https://www.clinisciences.com>



## Norway

Company: CliniSciences ApS  
Address: Oesterbrogade 226, st. 1,  
Copenhagen, 2100 - Denmark  
Telephone: +45 89 888 349  
Fax: +45 89 884 064  
Email: [norge@clinisciences.com](mailto:norge@clinisciences.com)  
Web: <https://www.clinisciences.com>



## Poland

Company: CliniSciences sp.Z.o.o.  
Address: ul. Rotmistrza Witolda Pileckiego 67  
lok. 200 - 02-781 Warszawa -Poland  
Telephone: +48 22 307 0535  
Fax: +48 22 307 0532  
Email: [polska@clinisciences.com](mailto:polska@clinisciences.com)  
Web: <https://www.clinisciences.com>



## Portugal

Company: Quimigen Unipessoal LDA  
Address: Rua Almada Negreiros, Lote 5, Loja 14,  
2615-275 Alverca Do Ribatejo - Portugal  
Telephone: +351 30 8808 050  
Fax: +351 30 8808 052  
Email: [info@quimigen.com](mailto:info@quimigen.com)  
Web: <https://www.quimigen.pt>



## Spain

Company: CliniSciences Lab Solutions  
Address: C/ Hermanos del Moral 13  
(Bajo E), 28019, Madrid - Spain  
Telephone: +34 91 269 40 65  
Fax: +34 91 269 40 74  
Email: [espana@clinisciences.com](mailto:espana@clinisciences.com)  
Web: <https://www.clinisciences.com>



## Sweden

Company: CliniSciences ApS  
Address: Oesterbrogade 226, st. 1,  
Copenhagen, 2100 - Denmark  
Telephone: +45 89 888 349  
Fax: +45 89 884 064  
Email: [sverige@clinisciences.com](mailto:sverige@clinisciences.com)  
Web: <https://www.clinisciences.com>



## Switzerland

Company: CliniSciences Limited  
Address: Marktgasse 18 8302 Kloten -  
Switzerland  
Telephone: +41 (044) 805 76 81  
Fax: +41 (044) 805 76 75  
Email: [switzerland@clinisciences.com](mailto:switzerland@clinisciences.com)  
Web: <https://www.clinisciences.com>



## UK

Company: CliniSciences Limited  
Address: 11 Progress Business center, Whittle  
Parkway, SL1 6DQ Slough- United Kingdom  
Telephone: +44 (0)1753 866 511  
or +44 (0) 330 684 0982  
Fax: +44 (0)1753 208 899  
Email: [uk@clinisciences.com](mailto:uk@clinisciences.com)  
Web: <https://www.clinisciences.com>



## USA

Company: Biotrend Chemicals LLC  
Address: c/o Carr Riggs Ingram,  
500 Grand Boulevard, Suite 210 Miramar  
Beach, FL 32550- USA  
Telephone: +1 850 650 7790  
Fax: +1 850 650 4383  
Email: [info@biotrend-usa.com](mailto:info@biotrend-usa.com)  
Web: <https://www.biotrend-usa.com>

