



- **NEW** improved heavy duty, kink-free EZBio® TPE V (Vacuum) tubing superior for performance under vacuum pressure
- **Vactrap2™** is available in autoclavable PP or chemically-resistant HDPE
- Each **Vactrap2™** includes the patented VersaCap closure for tangle-free opening and closing
- A 0.2µm PTFE inline vent filter retains any aerosolized bacteria and protects the house vacuum or pump
- The entire assembly is nested neatly in a secondary container with a small footprint and incredible stability against tipping
- Pinch clamp allows for vacuum adjustment
- 40 CFR 264.175 compliant
- Rated for continuous vacuum up to 1 Bar

Caution

- Do not use bottles at temperature greater than the rated maximum
- Do not flame bottles
- Do not use for long-term storage of strong organic compounds

Designed in USA



ISO 13485

The Vactrap2™ is a registered trademark of Foxx Life Sciences.


Foxx Life Sciences' **Vactrap2™ Round** is the world's most advanced vacuum trap system. Ideal for the protection of vacuum lines and scientists alike, **Vactrap2™ Round** safely contains waste from supernatant removal, chemical separation, and cell or tissue culture media aspiration. Best of all, **Vactrap2™ Round** comes ready to use; simply attach to your vacuum and begin your work.

The full **Vactrap2™ Round** system includes Foxx Heavy Duty Vacuum Bottle (autoclavable PP), each with a patented VersaCap® closure for tangle-free opening, no need to disconnect tubing, ever. NEW Improved heavy duty EZBio® TPE V tubing to avoid kinking. A 0.2µm PTFE inline vent filter retains any aerosolized bacteria, and the entire assembly is nestled neatly in a secondary container with a small footprint and incredible stability. A pinch clamp allows for the modulation of vacuum. Other **Vactrap2™ Round** options offer one bottle or setups without secondary containment. Spare EZFlow® Heavy Duty tubing with a 50mm vent filter are also available.

- 1 Strong, chemically-resistant, EZBio® TPE V tubing accepts Pasteur pipettes and other attachments for versatile aspiration
- 2 Patented VersaCap® technology eliminates tube tangling
- 3 Thick-walled, shatterproof heavy duty bottles are available in HDPE and PP options
- 4 Pinch clamp allows for modulation of vacuum at point of aspiration
- 5 Inline 0.2µm PTFE vent filter retains all aerosolized bacteria
- 6 Four feet of EZBio® TPE V tubing for flexibility in BSC or lab bench
- 7 Sturdy secondary container prevents spills and substantially decreases risk of tipping



Distributed by:
CliniSciences Group

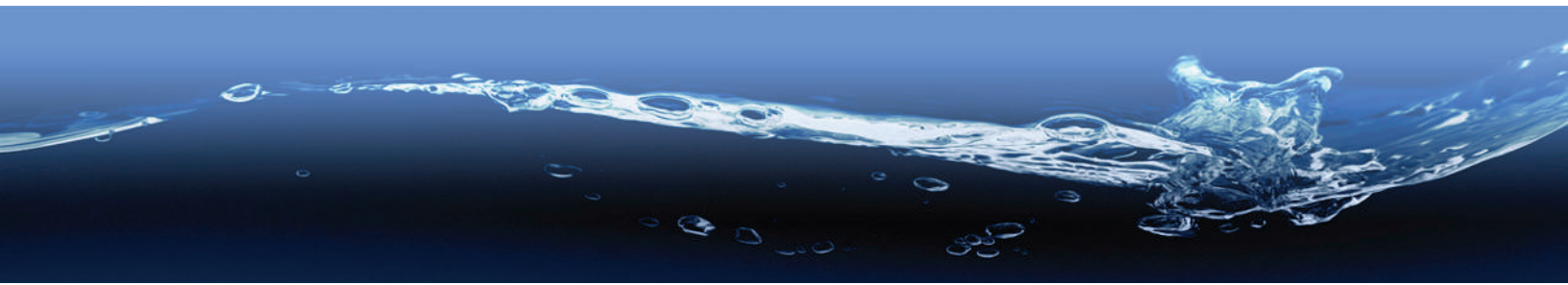

Part Number	Description	Material	Type of Assembly	U/M
305-1003-FLS	Vacrap2™ Round, 1L, PP, Red Bin, 1/4" ID Tubing	PP	Full Assembly	1/EA
305-1103-FLS	Vacrap2™ Round, 1L + 1L, PP, Red Bin, 1/4" ID Tubing	PP	Full Assembly	1/EA
305-2003-FLS	Vacrap2™ Round, 2L, PP, Red Bin, 1/4" ID Tubing	PP	Full Assembly	1/EA
305-2103-FLS	Vacrap2™ Round, 2 L + 1 L, PP, Red Bin, 1/4" ID Tubing	PP	Full Assembly	1/EA
305-4003-FLS	Vacrap2™ Round, 4 L, PP, Red Bin, 1/4" ID Tubing	PP	Full Assembly	1/EA
305-4103-FLS	Vacrap™ Round XL, 4L + 1L, PP, Red Bin, 1/4" ID Tubing	PP	Full Assembly	1/EA
305-4203-FLS	Vacrap™ Round XL, 4L + 2L, PP, Red Bin, 1/4" ID Tubing	PP	Full Assembly	1/EA
305-4403-FLS	Vacrap™ Round XL, 4L + 4L, PP, Red Bin, 1/4" ID Tubing	PP	Full Assembly	1/EA
305-2004-FLS	Vacrap2™ Round, 2L, PP, No Bin, 1/4" ID Tubing	PP	No Bin	1/EA
305-4004-FLS	Vacrap2™ Round, 4 L, PP, No Bin, 1/4" ID Tubing	PP	No Bin	1/EA

CliniSciences Group

Austria

Company: CliniSciences GmbH
Address: Sternwartestrasse 76, A-1180
Wien - Austria
Telephone: +43 720 115 580
Fax: +43 720 115 577
Email: oessterreich@clinisciences.com
Web: <https://www.clinisciences.com>



Belgium

Company: CliniSciences S.R.L
Address: Avenue Stalingrad 52, 1000
Brussels - Belgium
Telephone: +32 2 31 50 800
Fax: +32 2 31 50 801
Email: belgium@clinisciences.com
Web: <https://www.clinisciences.com>



Denmark

Company: CliniSciences ApS
Address: Oesterbrogade 226, st. 1,
Copenhagen, 2100 - Denmark
Telephone: +45 89 888 349
Fax: +45 89 884 064
Email: denmark@clinisciences.com
Web: <https://www.clinisciences.com>



Finland

Company: CliniSciences ApS
Address: Oesterbrogade 226, st. 1,
Copenhagen, 2100 - Denmark
Telephone: +45 89 888 349
Fax: +45 89 884 064
Email: suomi@clinisciences.com
Web: <https://www.clinisciences.com>



France

Company: CliniSciences S.A.S
Address: 74 Rue des Suisses, 92000
Nanterre- France
Telephone: +33 9 77 40 09 09
Fax: +33 9 77 40 10 11
Email: info@clinisciences.com
Web: <https://www.clinisciences.com>



Germany

Company: Biotrend Chemikalien GmbH
Address: Wilhelm-Mauser-Str. 41-43,
50827 Köln - Germany
Telephone: +49 221 9498 320
Fax: +49 221 9498 325
Email: info@biotrend.com
Web: <https://www.biotrend.com>



Iceland

Company: CliniSciences ApS
Address: Oesterbrogade 226, st. 1,
Copenhagen, 2100 - Denmark
Telephone: +45 89 888 349
Fax: +45 89 884 064
Email: island@clinisciences.com
Web: <https://www.clinisciences.com>



Ireland

Company: CliniSciences Limited
Address: Ground Floor, 71 lower Baggot street
Dublin D02 P593 - Ireland
Telephone: +353 1 6971 146
Fax: +353 1 6971 147
Email: ireland@clinisciences.com
Web: <https://www.clinisciences.com>



Italy

Company: CliniSciences S.r.l
Address: Via Maremmana inferiore 378
Roma 00012 Guidonia Montecelio - Italy
Telephone: +39 06 94 80 56 71
Fax: +39 06 94 80 00 21
Email: italia@clinisciences.com
Web: <https://www.clinisciences.com>



Netherlands

Company: CliniSciences B.V.
Address: Kraijenhoffstraat 137A,
1018RG Amsterdam, - Netherlands
Telephone: +31 85 2082 351
Fax: +31 85 2082 353
Email: nederland@clinisciences.com
Web: <https://www.clinisciences.com>



Norway

Company: CliniSciences AS
Address: c/o MerVerdi Munkerudtunet 10
1164 Oslo - Norway
Telephone: +47 21 988 882
Email: norge@clinisciences.com
Web: <https://www.clinisciences.com>



Poland

Company: CliniSciences sp.Z.o.o.
Address: ul. Rotmistrza Witolda Pileckiego 67
lok. 200 - 02-781 Warszawa -Poland
Telephone: +48 22 307 0535
Fax: +48 22 307 0532
Email: polaska@clinisciences.com
Web: <https://www.clinisciences.com>



Portugal

Company: Quimigen Unipessoal LDA
Address: Rua Almada Negreiros, Lote 5, Loja 14,
2615-275 Alverca Do Ribatejo - Portugal
Telephone: +351 30 8808 050
Fax: +351 30 8808 052
Email: info@quimigen.com
Web: <https://www.quimigen.pt>



Spain

Company: CliniSciences Lab Solutions
Address: C/ Hermanos del Moral 13
(Bajo E), 28019, Madrid - Spain
Telephone: +34 916 750 700
Fax: +34 91 269 40 74
Email: espana@clinisciences.com
Web: <https://www.clinisciences.com>



Sweden

Company: CliniSciences ApS
Address: Oesterbrogade 226, st. 1,
Copenhagen, 2100 - Denmark
Telephone: +45 89 888 349
Fax: +45 89 884 064
Email: sverige@clinisciences.com
Web: <https://www.clinisciences.com>



Switzerland

Company: CliniSciences AG
Address: Fracht Ost Flughafen Kloten
CH-8058 Zürich - Switzerland
Telephone: +41 (044) 805 76 81
Fax: +41 (044) 805 76 75
Email: switzerland@clinisciences.com
Web: <https://www.clinisciences.com>



UK

Company: CliniSciences Limited
Address: 11 Progress Business center, Whittle
Parkway, SL1 6DQ Slough- United Kingdom
Telephone: +44 (0)1753 866 511
or +44 (0) 330 684 0982
Fax: +44 (0)1753 208 899
Email: uk@clinisciences.com
Web: <https://www.clinisciences.com>



USA

Company: Biotrend Chemicals LLC
Address: c/o Carr Riggs Ingram,
500 Grand Boulevard, Suite 210 Miramar
Beach, FL 32550- USA
Telephone: +1 850 650 7790
Fax: +1 850 650 4383
Email: info@biotrend-usa.com
Web: <https://www.biotrend-usa.com>

